

**F&P** Eson™ 2

# NASAL MASK



## Confidence comes easy with Eson 2

F&P Eson 2 was designed to meet the needs of patients and sleep professionals, at every important milestone in the CPAP therapy journey.

OVER 20 DESIGN  
IMPROVEMENTS  
COMPARED TO  
F&P ESON!



# Confidence comes easy... at every CPAP therapy milestone

To improve on a mask designed for performance and comfort, our development team focused on ease of use, providing you and your patient confidence at every step from titration through to adherence.

## CPAP therapy journey

### FOR THE SLEEP PROFESSIONAL



1

#### SLEEP TEST

Selecting a nasal mask for titration

F&P Eson 2 retains the same RollFit™ Seal technology as F&P Eson, designed to deliver an effective seal to lessen the likelihood of leaks.



2

#### FIRST PATIENT CONTACT

Selecting a patient's mask for home use

The Intuitive headgear provides fitting options for patients, including the ability to stretch the mask directly over the head.



3

#### SET UP

Fitting and demonstrating the mask

Blue highlights – VisiBlue™ – incorporated into the mask provide visual cues for simpler, quicker and more intuitive mask education.



4

#### USE

Patient begins using mask

A design focus on ease of fitting and ease of use lessens the likelihood of call-backs and improves patient confidence.



5

#### ADHERENCE


Long-term adherence

Eson 2 has been designed to not just work on the first night, but to increase the chance of long-term therapy success.



Eson™ 2

# Key features and benefits

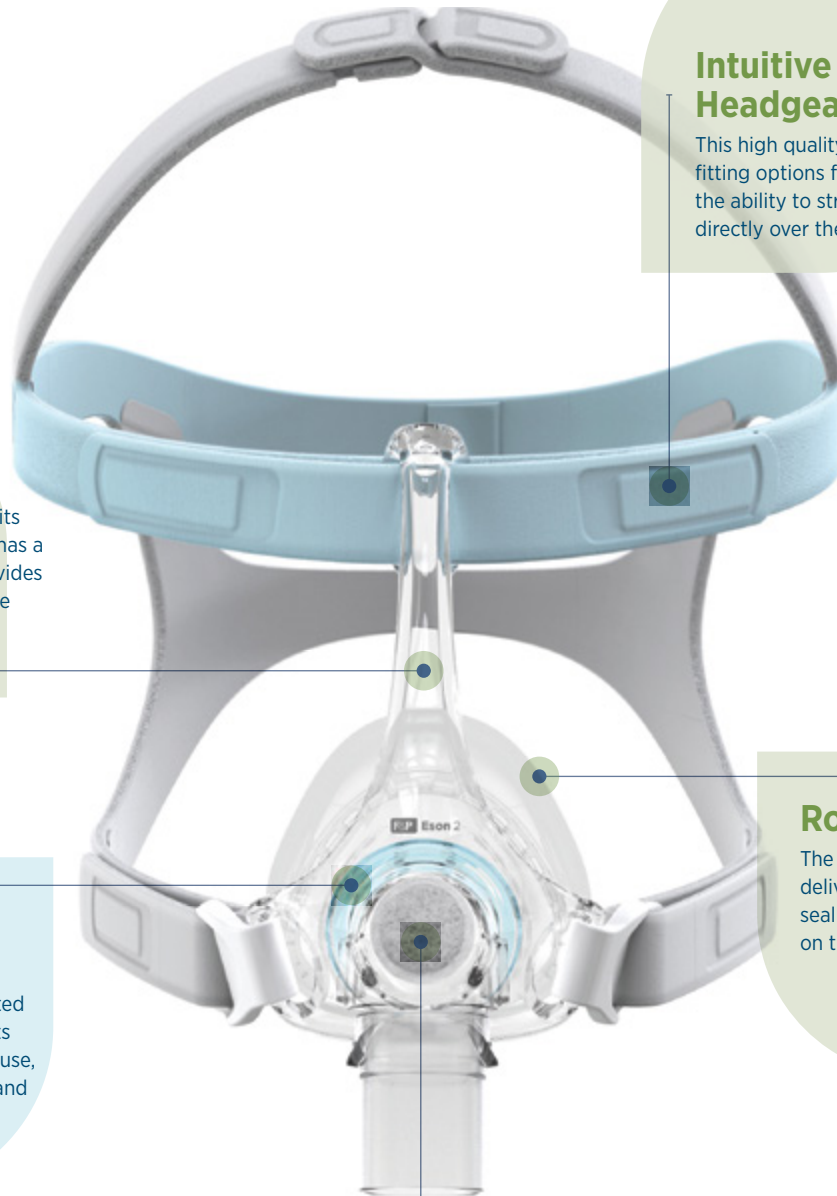

## Intuitive Headgear

This high quality headgear provides fitting options for patients, including the ability to stretch the mask directly over the head.

## Easy Frame

The Easy Frame - which fits all three seal sizes - now has a narrower profile. This provides a clearer line of sight while retaining seal stability.



## VisiBlue™

Blue highlights incorporated into key mask components assist patients with mask use, orientation, disassembly and reassembly.

## RollFit Seal




The RollFit Seal is designed to deliver a precise fit and effective seal by automatically adjusting on the bridge of the nose.

## Washable Diffuser

The innovative diffuser, which reduces noise and draft, is now washable and more durable.

# Comparing F&P Eson and F&P Eson 2

There have been over 20 design improvements made to make Fisher & Paykel Healthcare's best performing nasal mask even better.

COMPONENT	FEATURE	F&P ESON	F&P ESON2
<b>RollFit Seal</b>			
 <p>The RollFit Seal “rolls” back and forth, relieving pressure on the bridge of the nose and allowing adjustment to a wide variety of face types. The seal now includes VisiBlue highlights to aid with correct use, orientation, disassembly and reassembly.</p>	Larger seal alignment key		●
	Lighter and more compact seal		●
	VisiBlue highlights		●
<b>Intuitive Headgear</b>			
 <p>The Intuitive Headgear – incorporating the blue top strap to aid with orientation – can be fit multiple ways including stretching over the head, by using the traditional clip, or the fastener.</p>	ErgoFit Headgear		●
	Intuitive Headgear with VisiBlue highlights		●
	Improved Headgear strap fastener with 'stop tabs'		●
	Stitch-free, quality headgear		●
<b>Easy Frame</b>			
 <p>The Easy Frame fits all three seal sizes and now has a narrower profile for a clearer line of sight. The diffuser, swivel and clips are redesigned for enhanced ease of use.</p>	Removable ball-and-socket elbow for easier multi-patient cleaning		●
	One piece, washable diffuser		●
	Improved clip ergonomics		●
	Redesigned swivel with VisiBlue highlights		●

## NORTH AMERICA/ASIA-PACIFIC REORDER CODES

MASKS		400ESN221	ESON2 Headgear Small
ESN2SA	ESON2 Small	400ESN222	ESON2 Headgear Medium-Large
ESN2MA	ESON2 Medium	400ESN231	ESON2 Mask No H/G Small
ESN2LA	ESON2 Large	400ESN232	ESON2 Mask No H/G Medium
		400ESN233	ESON2 Mask No H/G Large
SPARE PARTS		400FPH241	F&P Swivel 2
400ESN211	ESON2 Seal Small	400ESN251	ESON2 Clips
400ESN212	ESON2 Seal Medium	400ESN261	ESON2 Diffuser
400ESN213	ESON2 Seal Large	400ESN271	ESON2 Elbow

## REST OF WORLD REORDER CODES

MASKS		400ESN221	ESON2 Headgear Small
ESN2SU	ESON2 Small	400ESN222	ESON2 Headgear Medium-Large
ESN2MU	ESON2 Medium	400FPH241	F&P Swivel 2
ESN2LU	ESON2 Large	400ESN251	ESON2 Clips
		400ESN261	ESON2 Diffuser
SPARE PARTS		400ESN271	ESON2 Elbow
400ESN211	ESON2 Seal Small		
400ESN212	ESON2 Seal Medium		
400ESN213	ESON2 Seal Large		